

RECONSTITUTING Gonal-f[®] RFF 75 IU VIAL

Prefilled syringe of Sterile Water for injection, USP (1 mL)



Gonal-f[®] RFF 75 IU Vial



27 G
1/2"

Injection Needle

18 G
1 1/2"


Mixing Needle



1. Flip the protective plastic cap off of the Gonal-f[®] RFF vial.



2. Wipe the top of the vial with an alcohol swab.

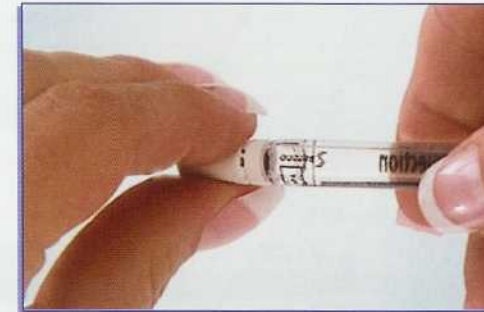

3. Hold the barrel of the prefilled syringe of Sterile Water in one hand. Firmly hold the plastic cap between the thumb and forefinger of the other hand and with a downward motion, gently snap and pull off the cap. If the gray cap remains, simply remove it.



4. Remove the safety seal cover of the 18 G 1 1/2" needle. Push the needle on the prefilled syringe until it is tightened. Holding the hub, or base, of the needle, secure the needle on the tip of the prefilled syringe and remove the needle cap.



5. With the vial of Gonal-f[®] RFF powder on a flat surface, insert the needle of the prefilled syringe straight down through the marked center circle of the rubber stopper. Slowly inject the water into the vial. DO NOT shake.



6. Invert the vial and pull back the 18 G 1 1/2" needle as far as needed and withdraw the entire contents of the vial. Remove the syringe from the vial.



7. If your dose requires more than one vial of Gonal-f[®] RFF 75 IU, use the mixture in the syringe to reconstitute the next vial of Gonal-f[®] RFF 75 IU powder. Use the same 18 G 1 1/2" needle and syringe to reconstitute additional vials.




8. Gently pull the plunger back to allow a small air space. Recap the needle. Twist and pull off needle from the syringe and discard in your sharps container.



9. Remove the safety seal cover of the 27 G 1/2" needle for injection. Push the needle on the prefilled syringe until it is tightened. Holding the hub, or base, of the needle, secure the needle on the tip of the prefilled syringe and remove the needle cap.



10. With the syringe pointing upward, gently tap on the syringe and slowly push the plunger until all air bubbles are gone and a drop of liquid appears on the tip of the needle.

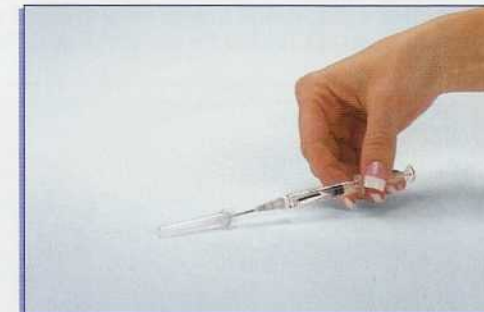

11. Recap the needle. The administration syringe is now filled with the prescribed dose of Gonal-f[®] RFF and is ready for administration. Use reconstituted Gonal-f[®] RFF 75 IU vial immediately.


RECONSTITUTING Gonal-f[®] MULTI-DOSE 450 IU VIAL




1. Flip the protective plastic cap off of the Gonal-f[®] Multi-Dose vial.



2. Wipe the top of the vial with an alcohol swab.

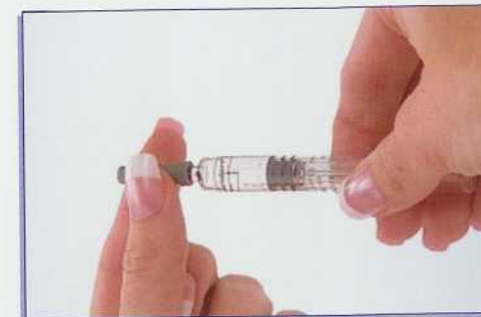

3. Carefully twist and pull off the rubber cap from the prefilled syringe of Bacteriostatic Water. To prevent contamination, avoid contact with the needle.

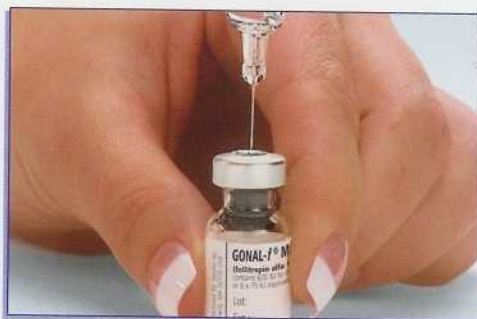

4. With the vial of Gonal-f[®] powder on a flat surface, insert the needle of the prefilled syringe straight down through the marked center circle of the rubber stopper. Slowly inject the water into the vial. DO NOT shake.



5. After all of the water has been injected into the vial, remove your finger from the plunger, allowing the plunger to rise to its original position. Withdraw the needle safely and dispose of it in a sharps container.



6. Wipe the top of the vial with an alcohol swab.



7. Remove the wrapper from the custom-dosing injection syringe. Carefully loosen and pull the plastic cap from the needle and avoid touching the needle.



8. With the vial of reconstituted Gonal-f[®] on a flat surface, insert the needle straight down through the marked center circle of the rubber stopper.



9. Without removing the needle from the vial, turn it upside down so that the needle points upward.




10. Slowly pull the plunger back until the syringe fills to slightly more than the unit marking that corresponds to your prescribed dose. Keeping the needle in the vial, slowly push the plunger to your prescribed dose. This will clear any air bubbles.



11. Carefully remove the syringe from the vial and recap the needle. The custom-dosing syringe is now filled with the prescribed dose of Gonal-f[®] and is ready for administration. Once reconstituted, store remaining solution at room temperature or refrigerate for up to 28 days.

SELF-INJECTING with Gonal-f® RFF Pen

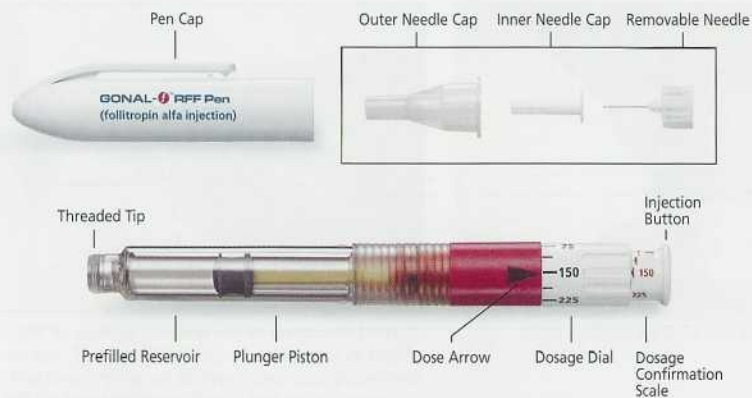

Priming the Pen: First use only



Follow instructions for Step 1: Attach needle. Set dose arrow at 37.5. Pull out injection button as far as it will go. Remove inner needle cap and hold pen so needle points upward. Tap prefilled reservoir to let any air bubbles rise to the top. (It's normal if a few small air bubbles remain.)

Keep needle pointing upward and push in injection button completely. Hold until you hear a click. A small amount of liquid should appear from needle tip. (If no liquid comes out, reprime the pen until it does.) Replace inner needle cap.

Disposal



Hold pen firmly by drug reservoir and replace outer needle cap. Grip outer needle cap firmly and unscrew pen from needle by turning pen counter-clockwise. Dispose of needle in safety container.


Storage

After use, replace pen cap and store pen at room temperature or in refrigerator for up to 28 days.

4 Easy Steps

step 1

Attach needle

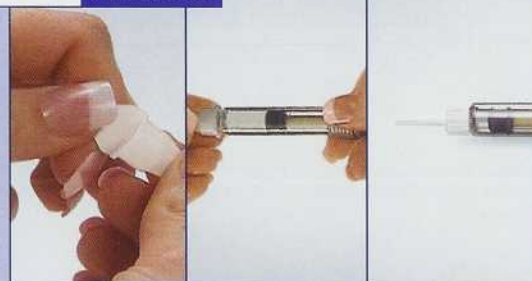

Remove pen cap and wipe threaded tip with alcohol wipe. Remove peel tab from outer needle cap. Press threaded tip of pen into open end of outer needle cap. Twist pen clockwise until needle is securely attached. Remove outer needle cap.

step 2

Set the dose




Turn dosage dial so your prescribed dose is lined up with the black dose arrow.

● **Recheck** that correct dose on dosage dial is lined up with black dose arrow before proceeding.

If the dialed dose is too high or too low, reset before pulling the injection button.

step 3

Load and confirm

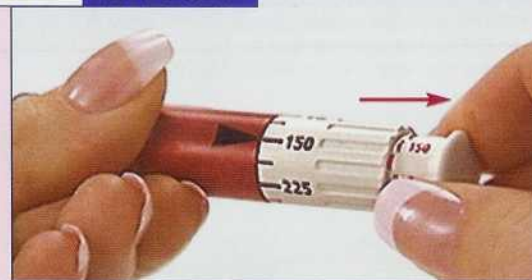

Load the dose by pulling out injection button as far as it will go. Confirm correct dose by looking at red dosage confirmation scale. The last visible flat red arrow indicates the loaded dose. (If loaded dose is lower than set dose, there is not enough drug to complete dose. Use new Gonal-f® RFF Pen to complete dose.)

NOTE: If set and loaded dose is too low: Turn dosage dial to correct dose and pull out injection button again.

If set and loaded dose is too high: Discard dose by pushing liquid out into safety container; repeat Step 2.

step 4

Inject



Clean recommended injection site with alcohol swab and let air dry. Remove inner needle cap. **Do not touch needle or allow it to touch any surface.** Inject needle into skin at a 90° angle and push injection button until it stops clicking. After last click, release pressure. Keep needle in the skin for at least 5 seconds. Remove needle and apply pressure as needed with gauze pad.